Draadloos massageapparaat

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Foreword

About this document

This manual contains all the information for correct, efficient and safe use of the product. Make sure you have fully read and understood the instructions in this manual before installing or using the product. Keep this manual for future use.

Original Instructions

This manual was originally written in Dutch. All other languages are translations of the original Dutch manual.

Service and warranty

If you need service or information regarding your product, please contact our customer service at service@auronic.nl.

We provide a 2-year warranty on our products. Customer must be able to show proof of purchase when claiming warranty. Product defects must be reported to us within 2 years from the date of purchase. Please note that the warranty does not apply: for damage caused by improper use or repair; for parts subject to wear and tear; for defects known to the customer at the time of purchase; for defects caused by the customer's own fault, for damage caused by third parties.

Table of contents

1. Introduction	14
1.1 Intended use	14
1.2 Product overview	14
1.3 Product specifications	15
1.4 Contents of package	15
2. Safety	16
3. Instructions for use	17
3.1 Operation	18
4. Troubleshooting	19
5. Maintenance and cleaning	19
6. Storage	19
7. Disposal and recycle	19
7.1 Disposal	20
7.2 Electric appliance	20
8. Disclaimer	20

1. Introduction

Thank you for choosing a Auronic product! Please always follow the instructions for safe usage.

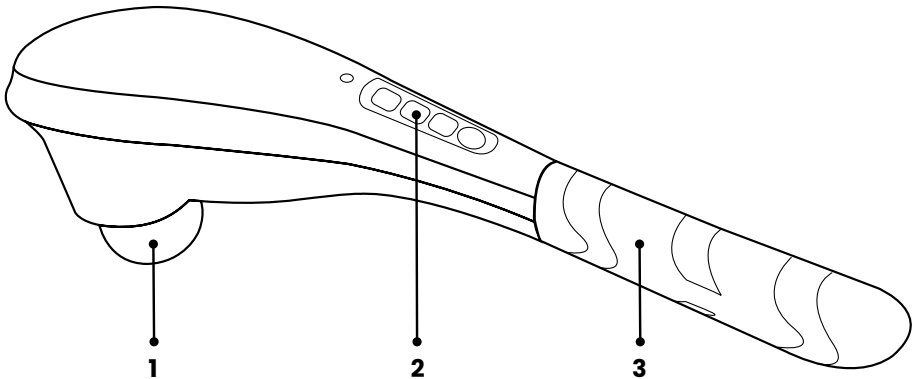
1.1 Intended use

A massage device is designed to relieve tension in muscles and improve blood flow. This device helps muscles recover after physical exertion. It also provides relief from pain by stimulating areas on the body that cause discomfort. The device can be used to reduce stress as it has a relaxing effect on the body and mind. It can be used for different parts of the body, such as back, legs, neck and shoulders. It contributes to well-being by promoting relaxation and increasing comfort.

The supplier is released from liability for any injury or damage caused by using the product beyond its intended use. Any modification of the product may affect your safety and warranty.

1.2 Product overview

1. Massage head
2. Control panel
3. Handle

1.3 Product specifications

Dimension	38,2 x 7,5 x 11,5 cm
Weight	1,15 kg
Voltage	220-240V AC
Current	0.5A
Frequency	50/60Hz
Input device	9V DC 1A
Power	12W
Battery	2 PCS 18650 2600mAh
Battery life	140 min when fully charged
Charging time	3.5 hours Please note that the massage device cannot be used while charging.
Pulses per minute	3600

1.4 Contents of package

When receiving your shipment, make sure that the contents are consistent with the packing list. Notify your distributor of any missing items. If the equipment appears to be damaged, file a claim immediately with the carrier and notify your distributor at once, giving a detailed description of any damage. Save the damaged parts and packing container to substantiate your claim.

- 1x Handheld massager
- 1x User manual
- 1x 3-point intensive massage head
- 1x 3-point massage head
- 1x Round massage head
- 1x curved massage head
- 1x 1-point massage head
- 1x charger

2. Safety

WARNING!

Please ensure that the packaging is kept out of the reach of children! Risk of suffocation!

Before use, read the manual and safety information carefully. Follow the instructions and keep the manual for future reference. The supplier is not responsible for failure to follow the safety information on the use and handling of the product or damage to property or personal injury.

DANGER

This category contains warnings for situations that could cause serious injury, damage or death if not avoided.

- Do not immerse the product, power cord or plug in water. If water does get inside or on the product, turn off the product immediately and contact our customer service department.
- To avoid electric shock, do not use the product if you suspect the power cord or plug is damaged. A damaged power cord or plug should be repaired or replaced by a professional (an authorized distributor).
- Do not unplug the product by pulling the power cord. This may cause fire or electric shock.
- If black smoke comes out of the product, unplug the product immediately. Wait until the smoke is gone before removing the contents of the product and contact our customer service department.
- If batteries/batteries leak, avoid contact with skin, eyes, open skin or wounds. Wear gloves. In case of accidental exposure to battery fluid, flush affected areas with plenty of clean water and seek medical attention immediately. Remove dead batteries/batteries from the product immediately as there is an increased risk of acid leakage.
- Do not break or disassemble the battery/battery pack.

**WARNING**

This category contains warnings for situations that could cause injury or damage if not avoided.

- Holding the product in the same place for too long can be painful and lead to burns.
- If the surface of the product is cracked, turn off the product to avoid electric shock.
- Fully charge the battery pack at least once a month to prevent battery damage.
- Do not charge the battery pack with any type of charger other than the charging station provided with the product.

**CAUTION**

This category contains warnings that require extra attention to avoid situations that could cause minor injuries or damage.

- Do not use the product in combination with massage oil, lotion, ointments, etc. This may cause damage to the product.

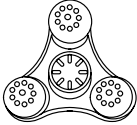
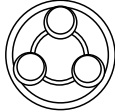
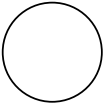
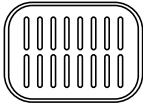
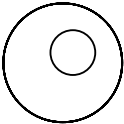
3. Instructions for use

**ATTENTION!**

- ◇ Charge the product before use. The charging light is red when the product is charging and green when it is charged.
- ◇ The product is equipped with an automatic emergency shutdown. After 20 minutes, the product automatically shuts off.
- ◇ If you apply too much pressure, the product automatically stops.

3.1 Operation

1. Select the desired massage head and attach it to the product.

	3-point intensive massage head Gentle, deep massage for a larger area.
	3-point massage head Gentle, deep massage.
	Round massage head Powerful general points massage.
	Curved massage head Arm and leg massage.
	1-point massage head Centered trigger-point massage.

2. Press the power button to turn the product on. The light will turn green when the product is on.
3. Select the desired mode by pressing the “mode” button and choose the desired speed.

Mode	Interval Mode Setting Setting 1: Continuous pulses Mode 2: 2 second pulses Setting 3: 1 second pulses Setting 4: 2 pulses per second
+	Increase massage force
-	Decrease massage force

4. Turn the product off if you are not using it (for an extended period of time). Charge the product every month to ensure longer battery life.

4. Troubleshooting

If the following problems occur while using the product, refer to the table below to find a solution. If the problems cannot be solved using these options, please contact our customer service.

Problem	Cause	Solution
The product does not work.	The battery is flat.	Charge the product. The power button lights green when the product is charged.
	The product has a malfunction.	Reset the product by turning it off, charging for about 3 seconds and turning it back on.
The product turns on, but does not massage.	It takes about 3 seconds for the product to start massaging.	Wait 3 seconds for the product to start massaging.
		If the product does not start massaging after 3 seconds, turn off the product and contact customer service.

5. Maintenance and cleaning

ATTENTION!

- ◇ Remove the product from the charger before cleaning it.

Clean the product with a damp, clean cloth and dry it thoroughly. Do not use abrasive or harsh cleaning agents.

6. Storage

Store the product in a clean and dry place. Do not store the product in extreme high or extreme low temperatures.

7. Disposal and recycle

Think about the environment and contribute to a cleaner living environment! You can contact your local council for all information on disposal options for discarded products.

7.1 Disposal

At the end of its life, do not dispose of the product with normal household waste, but hand it in at a government designated collection point for recycling.

Dispose of the packaging and accessories in accordance with the environmental regulations in force in your local council. Pay attention to the symbols on the different packaging materials and - if necessary - dispose of them separately. The packaging materials are marked with letters and numbers: 1-7: plastic, 20-22: paper and cardboard, 80-98 composite materials.

7.2 Electric appliance

The symbol shown here with the crossed-out garbage can indicates that this product is subject to European Directive (EU) 2012/19/EU. This directive states that at the end of its service life, you must not dispose of this product with normal household waste, but hand it in at a government designated collection point for recycling.

8. Disclaimer

Changes reserved; specifications are subject to change without stating grounds for doing so.

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