auronic

Waterflosser

NL/BE – Handleiding - Waterflosser

EN – User manual - Water flosser

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Foreword

About this document

This manual contains all the information for correct, efficient and safe use of the product. Make sure you have fully read and understood the instructions in this manual before installing or using the product. Keep this manual for future use.

Original Instructions

This manual was originally written in Dutch. All other languages are translations of the original Dutch manual.

Service and warranty

If you need service or information regarding your product, please contact our customer service at service@auronic.nl.

We provide a 2-year warranty on our products. Customer must be able to show proof of purchase when claiming warranty. Product defects must be reported to us within 2 years from the date of purchase. Please note that the warranty does not apply: for damage caused by improper use or repair; for parts subject to wear and tear; for defects

known to the customer at the time of purchase; for defects caused by the customer's own fault, for damage caused by third parties.

Table of contents

1. Introduction	24
1.1 Intended use	24
1.2 Product overview	25
1.3 Product specifications	26
1.4 Other features	26
1.5 Contents of package	26
2. Safety	27
2.1 General	27
2.2 Persons	29
2.3 Use	29
2.4 Medical use	30
2.5 Electronics	30
2.6 Batteries	31
3. Instructions for use	33
3.1 Before first use	33
3.2 Attach/remove attachment	33
3.3 Operation	34
3.4 After use	36
4. Maintenance and cleaning	37
5. Storage	37
6. Disposal and recycle	37
6.1 Disposal	37
6.2 Electric appliance	38
7. Symbols	39
8. Disclaimer	39

1. Introduction

Thank you for choosing a Auronic product! Please always follow the instructions for safe usage.


1.1 Intended use

A water flosser is a device that sprays pressurized water to remove plaque and food debris between teeth and along the gum line. The goal is to improve oral health by reaching areas that are difficult for toothbrushes and traditional floss to reach. This can help prevent gum disease and cavities. The device is especially useful for people with braces, implants or bridges, where plaque can easily accumulate. Regular use of a water flosser can help maintain healthy gums and clean teeth.

The supplier is released from liability for any injury or damage caused by using the product beyond its intended use. Any modification of the product may affect your safety and warranty.

1.2 Product overview

1. Attachment
2. Release button
3. Power Button
4. Mode button
5. Indicator light Normal
6. Indicator light Soft mode
7. Indicator light for massage mode (Pulse)
8. Battery indicator light
9. Water tank
10. Attachments
11. Charging point

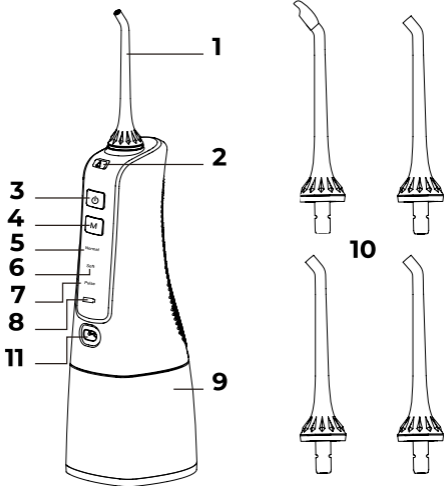1.3 Product specifications

Dimension	7.6x6.7x20.4cm
Weight	305g
Cord length	80cm
Voltage	5V DC
Current	1A
Battery	3.7V Lithium 1400mAh
Pressure	30-100psi
Battery Life	20 days (when fully charged)
IP value	IPX7
Pulses	1400-1800 times per minute
Water tank capacity	300ml
Charging time	4 hours
Battery life	20 days (if fully charged)
Modes	3: normal, soft and massage

1.4 Other features

- Removable water tank;
- 360° rotatable attachment;
- Wireless use;
- Memory function to remember your preset mode;
- 4 attachments including 1 suitable for bracket holders.

1.5 Contents of package

When receiving your shipment, make sure that the contents are consistent with the packing list. Notify your distributor of any missing items. If the equipment appears

to be damaged, file a claim immediately with the carrier and notify your distributor at once, giving a detailed description of any damage. Save the damaged parts and packing container to substantiate your claim.

- 1x Water flosser
- 1x User manual
- 4x Owner
- 1x Travel bag

2. Safety

WARNING!

Please ensure that the packaging is kept out of the reach of children! Risk of suffocation!

Before use, read the manual and safety information carefully. The supplier is not responsible for failure to follow the safety information on the use and handling of the product or damage to property or personal injury.

Always observe the following precautions when using 'the Water flosser', hereinafter referred to as 'the product'.

2.1 General

- If this is your first time using the product, it is important to read the manual and instructions carefully.
- If the product is malfunctioning, dropped or damaged, do not use the product and consult our customer service.

- Before each use, check the parts and accessories for any damage or wear (cracks/breakage) that may adversely affect performance, and if so, replace them with original or manufacturer-recommended parts and accessories.
- Do not use accessories not offered by the supplier. They may pose a danger to the user and damage the product.
- Use the product only for the purposes indicated in this manual.
- Do not disassemble, repair or restore the product. Doing so may cause damage to the product. For your own safety, these parts should only be replaced by authorized service personnel. Always send the product to a service center for repairs to prevent damage and personal injury.
- Be aware of the risks and consequences associated with improper use of the product. Misuse of the product may result in personal injury and damage to the product.
- Before first use, remove all packaging material and any stickers.
- The product is designed and intended for private use. Do not use the product for commercial purposes.
- Keep the product away from heat, areas where the temperature is relatively high, and heat sources such as ovens.
- Do not leave the product unattended when in use. When you are not using the product, turn it off and unplug it.
- Do not cover the product even when it is off.
- Fill the product with clean water. Do not put any other liquids in the product.

2.2 Persons

- This product may be used by children 8 years of age and older, provided that they are properly instructed and informed as to the safe use of the product and understand its potential hazards. Cleaning and maintenance of the product should not be performed by children unless they are over 8 years of age and supervised. Keep the product and cord out of the reach of children under 8 years of age.
- The product can be used by people with physical or mental disabilities, provided they are properly supervised and instructed in the safe use of the product and understand its potential hazards.
- Cleaning and maintenance of the product should not be performed by children unless they are over 8 years old and supervised.
- Do not allow children to play with the product. Use and store the product out of the reach of children.

2.3 Use

- Do not cover the product when it is on. Never place pillows or blankets on or against the product. This may cause a fire hazard, electric shock, personal injury or damage to the product.
- Do not clean the product while it is still plugged in.
- Do not use the product for more than 2 minutes at a time.
- Do not use the product to spray water into your ears, nose, or other sensitive areas. The water pressure may cause personal injury in these areas.
- Do not use the product near oxygen or aerosol sprays.
- Do not use the product when wearing oral jewelry.

2.4 Medical use

- If you experience pain or discomfort while using the product, stop immediately and consult your physician.
- This product is not a substitute for treatments from a doctor or physical therapist. Chronic symptoms may worsen.
- Do not use the product if you have a wound or inflammation on your tongue or in your mouth.
- Consult your doctor before use if you have dental conditions or have had oral surgery.
- This product does not serve as a substitute for treatments from a (dental) doctor.

2.5 Electronics

- If the power fails during use, turn off the product immediately and unplug the product.
- If the surface of the product is cracked, turn off the product to avoid electric shock.
- Check both power cord and plug regularly for signs of damage or wear, and if there is any, do not use the product and have the power cord and/or plug replaced or repaired by authorized service personnel.
- Use only the adapter and/or accessories provided by the supplier to prevent overheating, explosions, electric shock, fire, other damage and personal injury.
- To avoid electric shock, do not use the product if you suspect the power cord or plug is damaged. A damaged power cord or plug should be repaired or replaced by a professional (an authorized distributor).
- Unplug and turn off the product when not in use and before cleaning the product.

- Do not unplug the product by pulling the power cord. This may cause fire or electric shock.
- Do not let the cable hang over the edge of a table or countertop.
- Do not drop the product or subject it to large shocks.
- Do not place the cable on hot surfaces or near oil.
- Do not tuck the power cable under rugs or blankets. Make sure the cable is never in a walkway.
- Do not move the product by pulling the cord or twist the cord.
- Before first use, remove all packaging material and any stickers.
- Before connecting the product, make sure that the mains voltage indicated on the rating plate matches that of your mains supply.
- Make sure you do not trip over the cable and that the cable is not kinked, twisted, or stuck.
- If black smoke comes out of the product, unplug it immediately. Wait until the smoke is gone before removing the contents of the product and contact our customer service department.
- Do not insert objects into the openings of the product. Make sure that the moving parts are not blocked and that nothing is trapped.

2.6 Batteries

- Do not expose batteries/batteries to high temperatures or excessive heat.
- If batteries/batteries leak, avoid contact with skin, eyes, open skin or wounds. Wear gloves. In case of accidental exposure to battery fluid, flush affected

areas with plenty of clean water and seek medical attention immediately. Remove dead batteries/batteries from the product immediately, as there is an increased risk of acid leakage.

- Batteries/batteries must be inserted according to the correct polarity as indicated on the product.
- Batteries/batteries are deadly if swallowed. Keep batteries/batteries away from children.
- The batteries/batteries pose a choking hazard, keep them away from children and pets.
- The service life of batteries/batteries is approximately one year. Replace batteries/batteries on time.
- Use only batteries/batteries of the same or equivalent type.
- Do not mix old and new batteries/batteries.
- Use gloves to remove leaking batteries/batteries. Clean the battery compartment with a dry cloth.
- Do not dispose of batteries/batteries with household waste. Batteries/batteries contain chemicals that are harmful to the environment.
- Do not throw batteries/batteries into a fire. Do not attempt to crush, puncture or otherwise disassemble batteries/batteries for any reason.
- Leaking batteries/batteries may damage the product. Remove the batteries/batteries when not using the product for an extended period of time.
- Improper use of the batteries/batteries may cause explosion and leakage of battery fluid.
- Do not attempt to short-circuit the terminals or batteries/batteries.
- Do not attempt to charge non-rechargeable batteries/batteries. Rechargeable batteries/batteries must be removed from the product before charging.

- Recycle batteries/batteries if necessary. Return batteries/batteries to a government approved collection point for recycling.

3. Instructions for use

3.1 Before first use


Fully charge the product before the first use. This takes about 4 hours. Do not use the product while it is charging.

3.2 Attach/remove attachment

- Attach: place the attachment at the top of the opening of the product and turn until you hear a click.
- Remove: press the release button and remove the attachment from the opening by gently pulling it upward.

3.3 Operation


1. Make sure the product is off. Open the lid of the water reservoir and pour in water. Note: Make sure not to exceed the 300 ml line.
2. Choose the desired mode by pressing the mode button: standard, soft or massage. You can choose the desired mode both when the product is on and when the product is off.
3. Holding the product at an angle downward, place the end of the attachment in your mouth and aim at the teeth. Position the product so that the water sprays at a 90° angle to your gums. You can rotate the attachment to find the right position. Close your mouth slightly to avoid water spraying out.
4. Turn on the product by pressing the on/off button.
5. Start on the outside of your upper teeth (1). Then continue with the inside (2). Then continue with the inside of your lower teeth (3). Finally, proceed with the outside of your lower teeth (4).
6. The product automatically turns off after 2 minutes of continuous use.

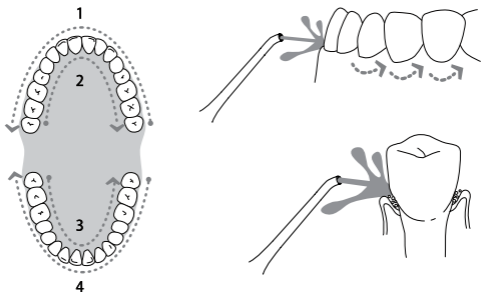- When the product is in massage mode, the battery life is displayed after 3 seconds.
- The battery is low if all three indicator lights flash during use. In this case, charge the product.
- The product has a memory function. When you turn on the product, it is automatically on the last used setting.

3.4 After use

1. Remove the water tank by holding the product with your left hand and the water tank with your right hand. Turn the water tank counterclockwise while gently pulling it down.
2. Empty the water tank by opening the water tank cover and pouring it empty.
3. Replace the water tank by turning it back clockwise on the product.

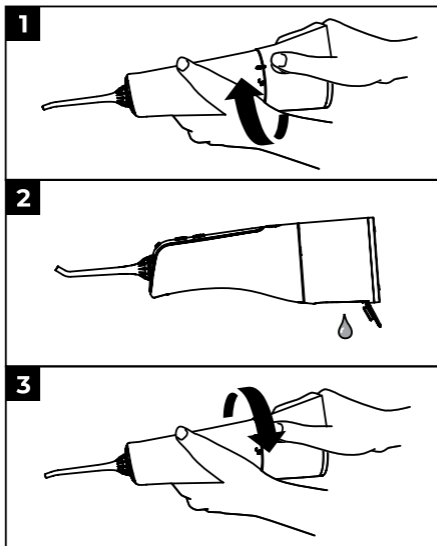4. Maintenance and cleaning

Clean the product with a damp, clean cloth and dry it well. Do not use abrasive or aggressive cleaning agents.

- If you have not used the product for a long time, let it fully charge before using it again. Charge the product in a timely manner to extend the battery life.

5. Storage

Store the product in a clean and dry place. Do not store the product in extreme high or extreme low temperatures.

6. Disposal and recycle

Think about the environment and contribute to a cleaner living environment! You can contact your local council for all information on disposal options for discarded products.

6.1 Disposal

At the end of its life, do not dispose of the product with normal household waste, but hand it in at a government designated collection point for recycling.

Dispose of the packaging and accessories in accordance with the environmental regulations in force in your local council. Pay attention to the symbols on the different packaging materials and - if necessary - dispose of them separately. The packaging materials are marked with letters and numbers: 1-7: plastic, 20-22: paper and cardboard, 80-98 composite materials.

6.2 Electric appliance

The symbol shown here with the crossed-out garbage can indicates that this product is subject to European Directive (EU) 2012/19/EU. This directive states that at the end of its service life, you must not dispose of this product with normal household waste, but hand it in at a government designated collection point for recycling.

7. Symbols

The following warnings and symbols are used in this manual, on the packaging and on the product (where applicable).

Symbol	Description
	Follow the instructions in the user manual and keep it for future reference.
	A warning with this symbol indicates a potentially hazardous situation that could result in serious damage and/or fatal injury.
	A note contains additional information relevant to the use of the product.
	A warning with this symbol indicates a potentially hazardous situation that may result in minor damage and/or moderate injury.
	Lithium battery/battery.
	Submersible: No water penetration if submerged (30 min at 1 m).

8. Disclaimer

Changes reserved; specifications are subject to change without stating grounds for doing so.

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©LifeGoods Group B.V.