# auronic

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---

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### FOREWORD

#### About this document

This manual contains all the information for correct, efficient and safe use of the product. Make sure you have fully read and understood the instructions in this manual before installing or using the product. Keep this manual for future use.

#### Original Instructions

This manual was originally written in Dutch. All other languages are translations of the original Dutch manual.

#### Service and warranty

If you need service or information regarding your product, please contact our customer service at [service@auronic.nl](mailto:service@auronic.nl).

We provide a 2-year warranty on our products. Customer must be able to show proof of purchase when claiming warranty. Product defects must be reported to us within 2 years from the date of purchase. Please note that the warranty does not apply: for damage caused by improper use or repair; for parts subject to wear and tear; for defects known to the customer at the time of purchase; for defects caused by the customer's own fault, for damage caused by third parties.

## TABLE OF CONTENTS

<b>1. Introduction</b>	<b>15</b>
1.1 Intended use	15
1.2 Product overview	15
1.3 Product specifications	16
1.4 Other features	16
1.5 Contents of package	16
1.6 Included accessories	16
<b>2. Safety</b>	<b>17</b>
2.1 General	17
2.2 People	18
2.3 Use	18
2.4 Medical remarks	18
2.5 Electrical	19
<b>3. Instructions for use</b>	<b>20</b>
3.1 Suitable for	20
3.2 Functionalities	20
<b>4. Maintenance and cleaning</b>	<b>21</b>
<b>5. Storage</b>	<b>21</b>
<b>6. Disposal and recycle</b>	<b>21</b>
6.1 Disposal	21
6.2 Electric appliance	21
<b>7. Symbols</b>	<b>22</b>
<b>8. CE Declaration of Conformity</b>	<b>23</b>
<b>9. Disclaimer</b>	<b>23</b>

## 1. INTRODUCTION

Thank you for choosing a Auronic product! Please always follow the instructions for safe usage.

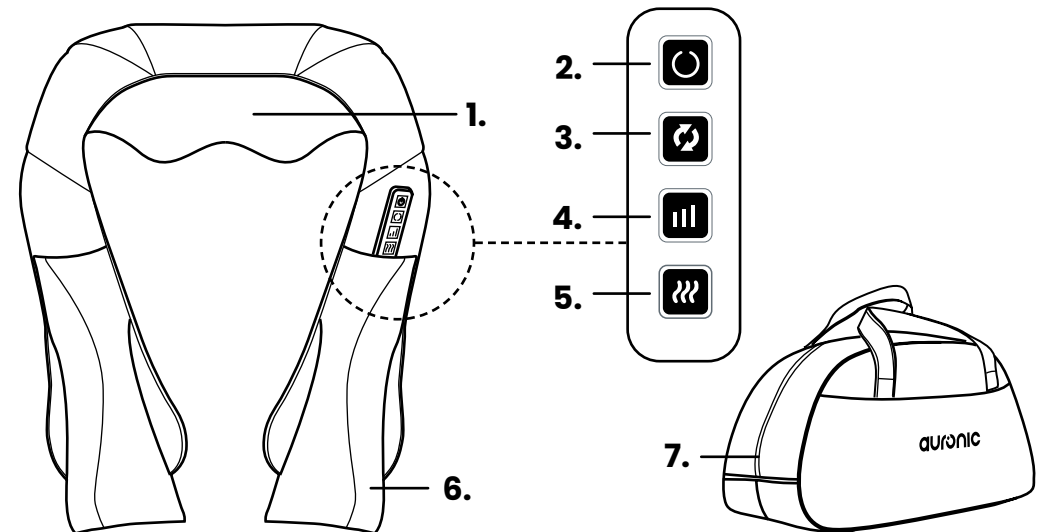
### 1.1 Intended use

A massager is designed to relieve tension in muscles and improve blood flow. This device helps muscles recover after physical exertion. It also provides relief from pain by stimulating areas on the body that cause discomfort. The device can be used to reduce stress as it has a relaxing effect on the body and mind. It can be used for different parts of the body, such as back, legs, neck and shoulders. It contributes to well-being by promoting relaxation and increasing comfort.

The supplier is released from liability for any injury or damage caused by using the product beyond its intended use. Any modification of the product may affect your safety and warranty.

### 1.2 Product overview

1. Massage heads
2. On/Off button
3. Selection button for the direction
4. Selection button for the speed
5. Selection knob for heat mode
6. Arm loop
7. Bag

**1.3 Product specifications**

Dimension	Product: 41 x 17 x 50 cm   Bag: 40 x 18 x 12 cm
Weight	1.4 kg
Color	AU3031: Black/Gray   AU3236: Black   AU3237: Black/Blue
Material	Product: Polyurethane (PU) leather   Bag: Leather
Safety class	0
Type of plug/adaptor	EU
Voltage	12V =
Amperage	5A
Frequency	50Hz
Power	24W
Overheating protection	Yes
Automatic switch off	Yes, after 15 min
Overheating protection	Yes
Massage type	Shiatsu
Heat function	Yes, up to 40°C
Number of massage heads	8
Number of speed settings	3

\*Recommended

**1.4 Other features**

- With adjustable massage direction.
- With arm loops to control the force yourself.
- With carrying bag to store the massager.

**1.5 Contents of package**

When receiving your shipment, make sure that the contents are consistent with the packing list. Notify your distributor of any missing items. If the equipment appears to be damaged, file a claim immediately with the carrier and notify your distributor at once, giving a detailed description of any damage. Save the damaged parts and packing container to substantiate your claim.

Auronic Massager with bag – Black/Gray – Model AU3032 – EAN 8720195254812

Auronic Massager with bag – Black – Model AU3249 – EAN 8720195259411

Auronic Massager with bag – Black/Blue – Model AU3250 – EAN 8720195259428

**1.6 Included accessories**

- 1x User manual
- 1x Power adapter
- 1x Bag

**2. SAFETY**

*Warning!*

*Please ensure that the packaging is kept out of the reach of children!  
Risk of suffocation!*

Before use, read the manual and safety information carefully. The supplier is not responsible for failure to follow the safety information on the use and handling of the product or damage to property or personal injury.

Always observe the following precautions when using ‘the Massager with bag’, hereinafter referred to as ‘the product’.

**2.1 General**

- If this is your first time using the product, it is important to read the user manual and instructions carefully.
- If the product is malfunctioning, dropped or damaged, do not use the product and consult our customer service.
- Before each use, check the parts and accessories for damage or wear (cracks/ breakage) that may adversely affect its operation and, if so, replace them with original or manufacturer-recommended parts and accessories.
- Do not use accessories not offered by the supplier. They may pose a danger to the user and damage the product.
- Use the product only for the purposes indicated in this user manual.
- Do not disassemble, repair or restore the product. Doing so may cause damage to the product. For your own safety, these parts should only be replaced by authorized service personnel. Always steer the product to a service center for repairs to prevent damage and personal injury.
- Be aware of the risks and consequences associated with improper use of the product. Misuse of the product may result in personal injury and damage to the product.
- Before first use, remove all packaging material and any stickers.
- The product is designed and intended for private use. Do not use the product for commercial purposes.
- Keep the product away from heat, areas where the temperature is relatively high and heat sources such as ovens.
- Do not insert objects into the openings of the product. Make sure that the moving parts are not blocked and that nothing is trapped.

**2.2 People**

- Keep the product and cord out of the reach of children.
- Do not let pets or small children chew on the power cord to avoid electric shock.
- Do not allow children to play with the product. Use and store the product out of the reach of children.
- The product should not be used by children under 8 years of age. Use of the product is allowed only when a supervisor is present and responsible for the user's safety.
- The product must not be used by people with physical or mental disabilities. Use of the product is permitted only in the presence of a supervisor responsible for the safety of the user.

**2.3 Use**

- Do not cover the product when it is on. Never place pillows or blankets on or against the product. This may cause a fire hazard, electric shock, personal injury or damage to the product.
- Parts of the product become hot. Do not clean the product while it is still hot. Wait until the product has cooled down before cleaning it.
- Do not use the product near a bathtub, shower, swimming pool, or near a (laundry) tub filled with water. Do not touch the product or power cord with wet hands. Do not allow the product to come into contact with water. Do not immerse the product in water. If water does get in or on the product, turn the product off immediately and contact our customer service department.
- Do not use the product for more than 15 minutes at a time.

**2.4 Medical remarks**

- Consult a doctor before using the product if you have unexplained pain, undergo medical treatment or use medical equipment.
- Do not use the product or consult a physician: if you are pregnant; if you have a pacemaker or other electronic implants; artificial joints, high blood pressure, heart disease; cerebral atherosclerosis, cerebral thrombosis, bone fractures, malignant tumors, circulation problems, varicose veins, open wounds, bruises, broken skin, inflammation of the veins, jaundice, diabetes, nervous disorders, acute inflammation or other health problems; during or immediately after major exertion; immediately after showering, taking a bath, or after the sauna; if you have been drinking alcohol; during or immediately after eating; if you have a fever; if you are in unexplained pain, undergoing medical treatment, or using medical equipment; if you have one or more of the following conditions: circulation problems, varicose veins, open wounds, bruises, broken skin, inflammation of the veins.
- If you experience pain or discomfort while using the product, stop immediately and consult your doctor.
- This user manual describes the possible uses of the product. Always seek the advice of a professional trainer or physician to adopt the correct posture for each different exercise. A professional trainer or doctor can reconcile your personal physical limits with the use of this product.

- This product is not a substitute for treatments from a physician or physical therapist. Chronic symptoms may worsen.

**2.5 Electrical**

- If the power fails during use, turn the product off immediately and unplug the product.
- If the surface of the product is cracked, turn off the product to avoid electric shock.
- Check both power cord and plug regularly for signs of damage or wear, and if there is any, do not use the product and have the power cord and/or plug replaced or repaired by authorized service personnel.
- Do not immerse the product, power cord or plug in water. If water does get inside or on the product, turn the product off immediately and contact our customer service department.
- Use and store the product so that the cables cannot be damaged.
- To avoid electric shock, do not use the product if you suspect the power cable or plug is damaged. A damaged power cord or plug should be repaired or replaced by a professional (an authorized distributor).
- Do not unplug the product by pulling the power cord. This may cause fire or electric shock.
- If black smoke comes out of the product, unplug the product immediately. Wait until the smoke is gone before removing the contents of the product and contact our customer service department.
- Do not leave the cable(s) hanging over the edge of a table or countertop.
- Do not place the cable(s) on hot surfaces or near oil.
- Do not tuck the cable(s) under rugs or blankets. Never lay the cable(s) in a walkway.
- Do not move the product by pulling the cord or twist the cord.
- Before connecting the product, make sure that the mains voltage indicated on the rating plate matches that of your mains supply.
- Make sure that you do not trip over the cable(s) and that the cable(s) are not kinked, twisted, or stuck.
- Do not leave the product unattended when in use. If you are not using the product for a long period of time, turn the product off and unplug it.

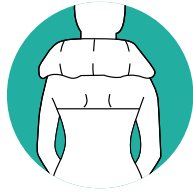
### 3. INSTRUCTIONS FOR USE

#### 3.1 Suitable for

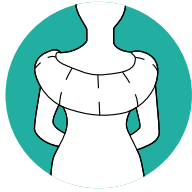
1. Neck and shoulders
2. Neck
3. Back
4. Middle
5. Abdomen
6. Upper legs
7. Legs
8. Feet
9. Calves

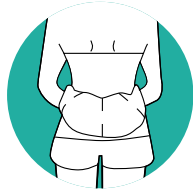
1.

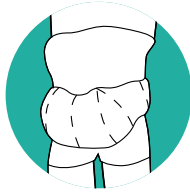
2.

3.

4.

5.

6.

7.

8.

9.

#### 3.2 Functionalities

	Place the product with the massage heads on the desired body part and start the massage by pressing the on/off button.
	Adjust the rotation direction of the massage heads by pressing the selection button. You can choose from counterclockwise or clockwise.
	Adjust the speed of the massage heads by pressing the selection button. You can choose from slow, medium and fast.
	Turn the heat mode of the massage heads on or off by pressing this button.

### 4. MAINTENANCE AND CLEANING

Unplug the power adapter cable from the product and power outlet before cleaning the product. Clean the product with a damp, clean cloth and dry it thoroughly. Do not use abrasive or harsh cleaning agents. Do not immerse the product in water or other liquids.

### 5. STORAGE

Store the product in a clean and dry place. Do not store the product in extreme high or extreme low temperatures.

### 6. DISPOSAL AND RECYCLE

Think about the environment and contribute to a cleaner living environment! You can contact your local council for all information on disposal options for discarded products.

#### 6.1 Disposal

At the end of its life, do not dispose of the product with normal household waste, but hand it in at a government designated collection point for recycling.

Dispose of the packaging and accessories in accordance with the environmental regulations in force in your local council. Pay attention to the symbols on the different packaging materials and - if necessary - dispose of them separately. The packaging materials are marked with letters and numbers: 1-7: plastic, 20-22: paper and cardboard, 80-98 composite materials.

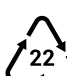
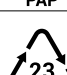
#### 6.2 Electric appliance

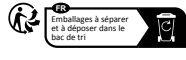
The symbol shown here with the crossed-out garbage can indicates that this product is subject to European Directive (EU) 2012/19/EU. This directive states that at the end of its service life, you must not dispose of this product with normal household waste, but hand it in at a government designated collection point for recycling.

## 7. SYMBOLS

The following warnings and symbols are used in this manual, on the packaging and on the product (where applicable).

Symbol	Description
	Follow the instructions in the user manual and keep it for future reference.
	A warning with this symbol indicates a potentially hazardous situation that could result in serious damage and/or fatal injury.
	A note contains additional information relevant to the use of the product.
	A warning with this symbol indicates a potentially hazardous situation that may result in minor damage and/or moderate injury.
	This product complies with European Directives 2014/30/EC (EMC - Electromagnetic Compatibility).
<b>LVD</b>	This product complies with 2014/35/EC (LVD - Low Voltage Directive) with regard to low voltage protection.
 <b>RoHS</b>	This product complies with all requirements and regulations in accordance with Directive 2011/65/EU better known as the RoHS (Restriction of Hazardous Substances) Directive. The directive contains regulations limiting the use of hazardous substances in electrical(on)ish equipment.
	Do not use the product if you have a pacemaker.
	Recycling symbol: Other plastics.
	Recycling symbol: Paper.
	Recycling symbol: Cardboard.

	Recycling symbol: This user manual can be recycled. (Only applicable in France)
	Recycling symbol: Separate packaging and put it in the sorting bin. (Only applicable in France)
	Recycling symbol: This product and its accessories can be recycled. (Only applicable in France)

## 8. CE DECLARATION OF CONFORMITY

This product complies with all requirements and regulations according to 2014/35/EU, Directive 2014/30/EU and Directive 2009/125/EC.

The complete EU Declaration of Conformity is available at the supplier.

## 9. DISCLAIMER

Changes reserved; specifications are subject to change without stating grounds for doing so.

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 LifeGoods B.V.