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NL/BE – Handleiding – Warmtekussen

EN – User manual – Heating Pad

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Foreword

About this document

This manual contains all the information for correct, efficient and safe use of the product. Make sure you have fully read and understood the instructions in this manual before installing or using the product. Keep this manual for future use.

Original Instructions

This manual was originally written in Dutch. All other languages are translations of the original Dutch manual.

Service and warranty

If you need service or information regarding your product, please contact our customer service at service@auronic.nl.

We provide a 2-year warranty on our products. Customer must be able to show proof of purchase when claiming warranty. Product defects must be reported to us within 2 years from the date of purchase. Please note that the warranty does not apply: for damage caused by improper use or repair; for parts subject to wear and tear; for defects known to the customer at the time of purchase; for defects caused by the customer's own fault, for damage caused by third parties.

Table of contents

1. Introduction	15
1.1 Intended use	15
1.2 Product overview	15
1.3 Product specifications	16
1.4 Other features	16
1.5 Contents of package	16
2. Safety	17
2.1 General	17
2.2 Persons	18
2.3 Use	18
2.4 Medical use	18
2.5 Electronics	18
3. Instructions for use	19
3.1 Operation	19
4. Maintenance and cleaning	19
5. Storage	20
6. Disposal and recycle	20
6.1 Disposal	20
6.2 Electric appliance	20
7. Symbols	21
8. CE Declaration of Conformity	22
9. Disclaimer	22

1. Introduction

Thank you for choosing a Auronic product! Please always follow the instructions for safe usage.


1.1 Intended use

An electric heat pad is designed to apply heat locally to specific parts of the body. It helps reduce stiffness and promote blood circulation. Because of its controlled heat delivery, it can help relax tense muscles. The pillow works through an electric element that produces heat, often adjustable to different heat levels.

The supplier is released from liability for any injury or damage caused by using the product beyond its intended use. Any modification of the product may affect your safety and warranty.

1.2 Product overview

1. Heat pad
2. Switch
3. Power cord

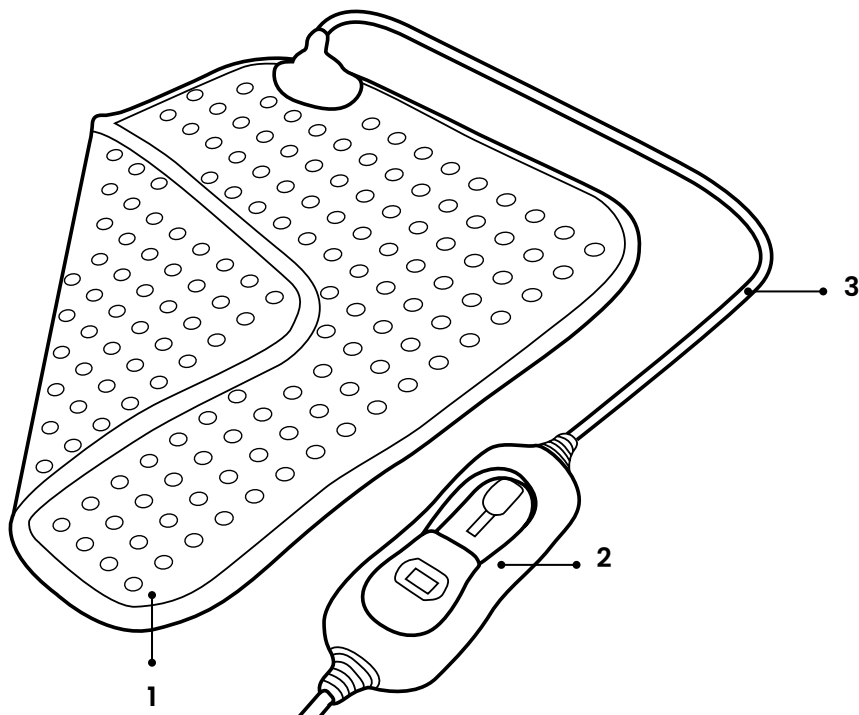1.3 Product specifications

Dimension	40 x 30 cm
Weight	670 gr
Cord length	240 cm
Material	Fleece
Voltage	220–240V AC
Frequency	50Hz
Power	100W
3 Heat settings	Low: +- 53 °C Medium: +- 59 °C High: +- 65 °C
Overheating protection	Yes
Automatic shut off	Yes, after 90 minutes of continuous use.

1.4 Other features

- Suitable for hand washing.
- Removable remote control with LED display.

1.5 Contents of package

When receiving your shipment, make sure that the contents are consistent with the packing list. Notify your distributor of any missing items. If the equipment appears to be damaged, file a claim immediately with the carrier and notify your distributor at once, giving a detailed description of any damage. Save the damaged parts and packing container to substantiate your claim.

- 1x Heating Pad
- 1x User manual

2. Safety

WARNING!

Please ensure that the packaging is kept out of the reach of children! Risk of suffocation!

Before use, read the manual and safety information carefully. The supplier is not responsible for failure to follow the safety information on the use and handling of the product or damage to property or personal injury.

Always observe the following precautions when using 'the Heating Pad', hereinafter referred to as 'the product'.

2.1 General

- If this is your first time using the product, it is important to read the manual and instructions carefully.
- If the product is malfunctioning, dropped or damaged, do not use the product and consult our customer service.
- Before each use, check the parts and accessories for damage or wear (cracks/ breakage) that may adversely affect its operation and, if so, replace them with original or manufacturer-recommended parts and accessories.
- Do not use accessories not offered by the supplier. They may pose a danger to the user and damage the product.
- Use the product only for the purposes indicated in this manual.
- Do not disassemble, repair or restore the product. Doing so may cause damage to the product. For your own safety, these parts should only be replaced by authorized service personnel. Always send the product to a service center for repairs to prevent damage and personal injury.
- Be aware of the risks and consequences associated with improper use of the product. Misuse of the product may result in personal injury and damage to the product.
- Before first use, remove all packaging material and any stickers.
- The product is suitable for indoor use only. Do not use the product outdoors. Doing so may result in damage to the product.
- Keep the product away from heat, areas where the temperature is relatively high, and heat sources such as ovens.
- Do not use the product if it is folded or creased.
- Keep the product away from open flames, lit cigarettes and other heat sources and/or hot surfaces.
- Do not expose the product to direct sunlight, extreme temperatures, humidity or wet environments for an extended period of time.
- Never place the product on or over the head.

2.2 Persons

- The product should not be used on persons who are insensitive to heat or who need a lot of protection because they cannot react to overheating (e.g., diabetics, persons with disease-induced skin changes or scars in the area of application, after taking analgesics or narcotics, or alcohol).
- Children under 3 years of age should not use this product because they cannot react to overheating.
- Keep the product and cord out of the reach of children.
- Do not let pets or small children chew on the power cord to avoid electric shock.
- Do not allow children to play with the product. Use and store the product out of the reach of children.

2.3 Use

- When using the product on an adjustable bed, make sure that the product cannot crease and the control panel cannot be covered when adjusting the bed.
- Do not use the product near objects or surfaces that are sensitive to heat.
- Parts of the product will become hot. Do not clean the product while it is still hot. Wait until the product has cooled down before cleaning it.

2.4 Medical use

- Do not use the product if you have circulatory problems and do not use the product on areas sensitive to heat. This can lead to burns!
- Do not use the product or consult a doctor before use if: you are pregnant; you have a pacemaker or other electronic implants; you have circulation problems, varicose veins, open wounds, bruises, broken skin, inflammation of the veins or other skin conditions; you are sensitive to heat.
- Do not use the product; on sunburned or sensitive skin; swollen or inflamed areas; areas with poor circulation, broken or injured skin, bruises or skin irritations; on unconscious or sleeping people; in the same place, in bed, for long periods of time. This may cause symptoms or pain.
- This product is not a substitute for treatments from a physician or physical therapist. Chronic symptoms may worsen.

2.5 Electronics

- Use only the supplied AC charger. Using any other charger may cause irreparable damage to the product and the environment and may be dangerous.
- If the power fails during use, immediately turn off the product and unplug it.
- If the surface of the product is cracked, turn off the product to avoid electric shock.
- Check both power cord and plug regularly for signs of damage or wear, and if there is any, do not use the product and have the power cord and/or plug replaced or repaired by authorized service personnel.
- Do not immerse the product, power cord or plug in water. If water does get inside or on the product, turn the product off immediately and contact our customer service department.

- Use and store the product so that the cables cannot be damaged.
- To avoid electric shock, do not use the product if you suspect the power cord or plug is damaged. A damaged power cord or plug should be repaired or replaced by a professional (an authorized distributor).
- Do not unplug the product by pulling the power cord. This may cause fire or electric shock.
- If black smoke comes out of the product, unplug the product immediately. Wait until the smoke is gone before removing the contents of the product and contact our customer service department.
- Do not leave the cable(s) hanging over the edge of a table or countertop.
- Do not place the cable(s) on hot surfaces or near oil.
- Do not tuck the cable(s) under rugs or blankets. Never lay the cable(s) in a walkway.
- Do not move the product by pulling the cord or twist the cord.
- Before connecting the product, make sure that the mains voltage indicated on the rating plate matches that of your mains supply.
- Make sure you do not trip over the cord(s) and that the cord(s) are not kinked, twisted, or stuck.

3. Instructions for use

ATTENTION!

- ◇ The product should only be used on the back and neck. Do not use the product near your eyes or other sensitive areas. Never place the product on or over the head. Do not sit on the product.
- ◇ Prolonged use on a high setting may cause burns.

3.1 Operation

1. Plug in and slide the switch from position 0 to 1. The product heats up slowly.
2. Slide the switch to position 2 or 3 for a higher temperature. Slide the switch back to a lower setting if the pillow gets too hot.
3. Slide the switch back to setting 1 before you go to sleep or if you use the product for a long time.
4. Flip the switch back to 0 to turn the product off.

4. Maintenance and cleaning

- Unplug the product and let it cool down for at least 10 minutes.
- Disconnect the plug and remote control from the product.
- Wash the product by hand.
- Use lukewarm water and mild detergent and squeeze the product gently.
- Rinse the product several times to rinse out all the detergent.
- Let the product air dry.
- Do not put the product in the dryer, use a hair dryer or other dryer, or iron the product.
- Use the product only when it is completely dry again.

5. Storage

Store the product in a clean and dry place. Do not store the product in extreme high or extreme low temperatures.

- Let the product cool completely before folding and storing it.
- Do not fold the product too tightly.
- Be careful not to damage the power cord when folding it.
- Do not place objects on the product when it is stored.

6. Disposal and recycle

Think about the environment and contribute to a cleaner living environment! You can contact your local council for all information on disposal options for discarded products.

6.1 Disposal

At the end of its life, do not dispose of the product with normal household waste, but hand it in at a government designated collection point for recycling.

Dispose of the packaging and accessories in accordance with the environmental regulations in force in your local council. Pay attention to the symbols on the different packaging materials and – if necessary – dispose of them separately. The packaging materials are marked with letters and numbers: 1-7: plastic, 20-22: paper and cardboard, 80-98 composite materials.

6.2 Electric appliance

The symbol shown here with the crossed-out garbage can indicates that this product is subject to European Directive (EU) 2012/19/EU. This directive states that at the end of its service life, you must not dispose of this product with normal household waste, but hand it in at a government designated collection point for recycling.

7. Symbols

The following warnings and symbols are used in this manual, on the packaging and on the product (where applicable).

Symbol	Description
	Follow the instructions in the user manual and keep it for future reference.
	A warning with this symbol indicates a potentially hazardous situation that could result in serious damage and/or fatal injury.
	A note contains additional information relevant to the use of the product.
	A warning with this symbol indicates a potentially hazardous situation that may result in minor damage and/or moderate injury.
	Suitable for indoor use.
	Wash by hand.
	Do not tumble dry.
	Do not bleach.
	Do not dry clean.
	Do not iron.
	Do not use if folded or rolled!
	Do not stitch with needles!

8. CE Declaration of Conformity

This product complies with all requirements and regulations according to 2014/35/EU, Directive 2014/30/EU and Directive 2009/125/EC.

The complete EU Declaration of Conformity is available at the supplier.

9. Disclaimer

Changes reserved; specifications are subject to change without stating grounds for doing so.

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©LifeGoods B.V.