Infrarood warmtekussen

NL/BE – Handleiding – Infrarood warmtekussen

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Foreword

About this document

This manual contains all the information for correct, efficient and safe use of the product. Make sure you have fully read and understood the instructions in this manual before installing or using the product. Keep this manual for future use.

Original Instructions

This manual was originally written in Dutch. All other languages are translations of the original Dutch manual.

Service and warranty

If you need service or information regarding your product, please contact our customer service at service@auronic.nl.

We provide a 2-year warranty on our products. Customer must be able to show proof of purchase when claiming warranty. Product defects must be reported to us within 2 years from the date of purchase. Please note that the warranty does not apply: for damage caused by improper use or repair; for parts subject to wear and tear; for defects known to the customer at the time of purchase; for defects caused by the customer's own fault, for damage caused by third parties.

Table of contents

1. Introduction	17
1.1 Intended use	17
1.2 Product overview	17
1.3 Product specifications	18
1.4 Material	18
1.5 Contents of package	19
2. Safety	19
2.1 General	19
2.2 Persons	20
2.3 Usage	20
2.4 Medical remarks	21
2.5 Electrical	21
2.6 Battery	22
3. Instructions for use	22
3.1 Charging	22
3.2 Battery	22
3.3 Heating up the product	23
3.4 Step-down function	23
3.5 Charging and using product	23
4. Maintenance and cleaning	23
5. Storage	23
6. Disposal and recycle	24
6.1 Disposal	24
6.2 Electric appliance	24
7. Symbols	25
8. CE Declaration of Conformity	25
9. Disclaimer	26

1. Introduction

Thank you for choosing a Auronic product! Please always follow the instructions for safe usage.


1.1 Intended use

An infrared heating pillow provides targeted heat and relaxation. It emits infrared radiation that penetrates deep into muscles and joints. This creates a pleasant, constant warmth that helps relieve tension and stress. The heating pillow is easy to use and suitable for use during moments of rest, such as on the couch or in bed. Thanks to the deep penetration of infrared radiation, the heating pillow provides effective heat treatment, promoting comfort and relaxation.

The supplier is released from liability for any injury or damage caused by using the product beyond its intended use. Any modification of the product may affect your safety and warranty.

1.2 Product overview

1. On/Off
2. Level indicator
3. Pillow

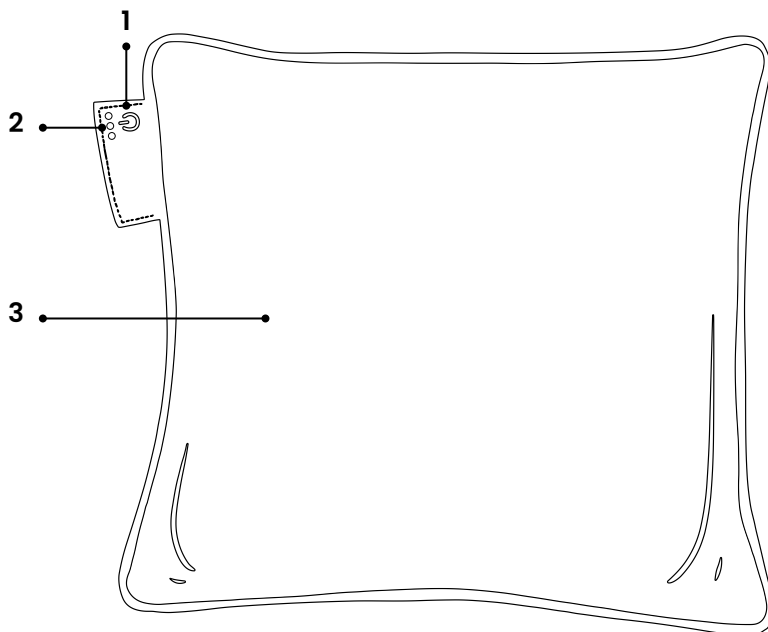1.3 Product specifications

	AU3337/AU3338/AU3339/ AU3340/AU3341/AU3342/ AU3343/AU3344/AU3345/ AU3346/AU3347/AU3348/ AU3349/AU3350/AU3351/ AU3352/AU3353/AU3354/ AU5052/AU5054/AU5056	AU3355/AU3356/AU3357/ AU3358/AU3359/AU3360/ AU3361/AU3362/AU3363/ AU3364/AU3365/AU3366/ AU3367/AU3368/AU3369/ AU3370/AU3371/AU3372/ AU5057/AU5059/AU5061
Dimension	45 x 45 cm	45 x 60 cm
Weight	650 g	680 g
Cord length	2.5 m	
Type of plug/adaptor	EU/AC	
Adapter voltage	Input: 100-240V AC Output: 8.4V DC 2.5A	
Frequency	50/60Hz	
Power	14.8W	
Power consumption	0.0148kWh	
Battery type	7.4V Lithium-ion	
Battery charging time*	±3 hours	
Battery capacity	5000mAh	
Replaceable battery	No	
Overheating protection	Yes	
Heating element	TPC fiber	

* This refers to the charging time for a full charge when the battery is depleted. Actual charging time may vary depending on ambient temperature and battery condition.

1.4 Material

- AU3337/AU3338/AU3339/AU3340/AU3341/AU3355/AU3356/AU3357/AU5052/AU5054/
AU5056/AU5057/AU5059/AU5061:
Side 1: 70% wool and 30% polyester / Side 2: 100% polyester
- AU3347/AU3348/AU3349/AU3350/AU3351/AU3365/AU3366/AU3367/AU3368/AU3369:
Side 1: 50% cotton and 50% polyester / Side 2: 100% polyester
- AU3342/AU3343/AU3344/AU3345/AU3346/AU3352/AU3353/AU3354/AU3358/AU3359/
AU3360/AU3361/AU3362/AU3363/AU3364/AU3370/AU3371/AU3372:
100% polyester

1.5 Contents of package

When receiving your shipment, make sure that the contents are consistent with the packing list. Notify your distributor of any missing items. If the equipment appears to be damaged, file a claim immediately with the carrier and notify your distributor at once, giving a detailed description of any damage. Save the damaged parts and packing container to substantiate your claim.

- 1x Infrared heating pillow
- 1x User manual
- 1x Cover
- 1x Battery
- 1x Charger

2. Safety

WARNING!

Please ensure that the packaging is kept out of the reach of children! Risk of suffocation!

Before use, read the manual and safety information carefully. The supplier is not responsible for failure to follow the safety information on the use and handling of the product or damage to property or personal injury.

Always observe the following precautions when using 'the Infrared heating pillow', hereinafter referred to as 'the product'.

2.1 General

- If this is your first time using the product, it is important to read the user manual and instructions carefully.
- If the product is malfunctioning, dropped or damaged, do not use the product and consult our customer service.
- Before each use, check the parts and accessories for damage or wear (cracks/ breakage) that may adversely affect its operation and, if so, replace them with original or manufacturer-recommended parts and accessories.
- Do not use accessories not offered by the supplier. They may pose a danger to the user and damage the product.
- Use the product only for the purposes indicated in this user manual.
- Do not disassemble, repair or restore the product. Doing so may cause damage to the product. For your own safety, these parts should only be replaced by authorized service personnel. Always steer the product to a service center for repairs to prevent damage and personal injury.

- Be aware of the risks and consequences associated with improper use of the product. Misuse of the product may result in personal injury and damage to the product.
- Before first use, remove all packaging material and any stickers.
- The product is suitable for indoor use only. Do not use the product outdoors. Doing so may result in damage to the product.
- Keep the product away from heat, areas where the temperature is relatively high, and heat sources such as ovens.
- Do not use the product if it is folded or creased.
- Keep the product away from open flames, lit cigarettes and other heat sources and/or hot surfaces.
- Do not expose the product to direct sunlight, extreme temperatures, humidity or wet environments for an extended period of time.
- Do not allow the product to get wet.
- Do not place other objects on the product when cleaning up and storing it.

2.2 Persons

- This product may be used by children 3 years of age and older provided they are properly supervised and instructed in the safe use of the product and understand its potential hazards. Cleaning and maintenance of the product should not be performed by children unless they are over 8 years of age and supervised.
- To avoid electric shock, do not allow pets or small children to chew on the power cord.
- Do not allow children to play with the product. Use and store the product out of the reach of children.
- Place the product out of reach of pets and other animals.
- The product should not be used on persons who are insensitive to heat or who need a lot of protection because they cannot react to overheating (e.g. diabetics, persons with disease-induced skin changes or scars in the area of application, after taking painkillers or narcotics or alcohol).
- Children under 3 years of age should not use this product because they cannot react to overheating.
- The product may be used by people with physical or mental disabilities, provided they are properly supervised and informed regarding the safe use of the product and understand its possible dangers.

2.3 Usage

- Do not use the product near a bathtub, shower, swimming pool, or near a (laundry) tub filled with water. Do not touch the product or power cord with wet hands. Do not allow the product to come into contact with water. Do not immerse the product in water. If water does get in or on the product, turn the product off immediately and contact our customer service department.
- If using the product on an adjustable bed, make sure the product cannot crease and the control panel cannot be covered when adjusting the bed.
- Do not use the product near objects or surfaces that are sensitive to heat.
- Keep the product away from extreme heat sources or direct sunlight.

- The product should not be used in places where special conditions prevail, such as the presence of a corrosive or explosive atmosphere.
- Do not sleep with the product, the user may suffer burns or heat stroke.

2.4 Medical remarks

- Do not use the product if you have circulatory problems and do not use the product in areas sensitive to heat. This may result in burns!
- Do not use the product or consult a doctor before use if: you are pregnant; you have a pacemaker or other electronic implants; you have circulation problems, varicose veins, open wounds, bruises, broken skin, inflammation of the veins or other skin conditions; you are sensitive to heat.
- The product is not intended for medical use in hospitals and/or health care facilities.

2.5 Electrical

- If the power fails during use, turn off the product immediately and unplug the product.
- If the surface of the product is cracked, turn off the product to avoid electric shock.
- Check both power cord and plug regularly for signs of damage or wear, and if there is any, do not use the product and have the power cord and/or plug replaced or repaired by authorized service personnel.
- Do not immerse the product, power cord or plug in water. If water does get inside or on the product, turn the product off immediately and contact our customer service department.
- Use and store the product so that the cables cannot be damaged.
- To avoid electric shock, do not use the product if you suspect the power cord or plug is damaged. A damaged power cord or plug should be repaired or replaced by a professional (an authorized distributor).
- Do not unplug the product by pulling the power cord. This may cause fire or electric shock.
- If black smoke comes out of the product, unplug the product immediately. Wait until the smoke is gone before removing the contents of the product and contact our customer service department.
- Do not leave the cable(s) hanging over the edge of a table or countertop.
- Do not place the cable(s) on hot surfaces or near oil.
- Do not tuck the cable(s) under rugs or blankets. Never lay the cable(s) in a walkway.
- Do not move the product by pulling the cord or twist the cord.
- Before connecting the product, make sure that the mains voltage indicated on the rating plate matches that of your mains supply.
- Make sure you do not trip over the cord(s) and that the cord(s) are not kinked, twisted, or stuck.
- Use only the AC charger provided. Use of any other charger may cause irreparable damage to the product and the environment and may be dangerous.
- Do not use an extension cord. Make sure the product is plugged into a grounded outlet to which no other devices are connected.
- The product should only be connected to a grounded 220-240V outlet.

2.6 Battery

- If the battery leaks, avoid contact with skin, eyes, open skin or wounds. Wear gloves. In case of accidental exposure to battery fluid, flush affected areas with plenty of clean water and seek medical attention immediately. Remove the empty battery pack from the product immediately, as there is an increased risk of acid leakage.
- Do not break or disassemble the battery pack.
- Always check the battery for signs of damage or leaks before using the product. Do not charge damaged or leaking batteries.
- The charger is designed for a specific voltage. Always check that the mains voltage matches the voltage on the rating plate.
- Battery life is best ensured by not charging the product until the battery is empty. Try to avoid intermediate charging times.
- Do not dispose of the battery into household waste. The battery contains chemicals that are harmful to the environment.
- Do not throw the battery into a fire. Do not attempt to crush, puncture or otherwise disassemble the battery pack for any reason.
- Remove the battery pack from the product first before disposing of the product.
- Fully charge the product before using it for the first time. This may take a few hours.
- Do not expose the battery pack to high temperatures or excessive heat.
- Do not heat the battery pack or place it near a heat source.

3. Instructions for use

3.1 Charging

1. Connect the charger to the back of the switch.
2. When the 3 indicators are lit continuously, the battery is fully charged.
3. When the product is not in use after fully charged, turn off the product and always unplug the power cord.

3.2 Battery

To check the battery status, press and hold the power button on the switch for 2 seconds. Over the next 5 seconds, 1 to 3 lights indicate the percentage of battery charge:

- 1 light: 0-30%
- 2 lights: 30-80%
- 3 lights: 80-100%

Based on the level selected, the battery will last for a specific quantity of hours:

- Level 1: 4 hours
- Level 2: 3 hours
- Level 3: 2.5 hours

**ATTENTION!**

- ◇ Only charge the original battery by connecting it to the product and only use the original charger to ensure safe and optimal performance of the product. Using any other battery may cause irreparable damage to the product and the environment and may be dangerous.
- ◇ Actual battery life may vary depending on ambient temperature and battery condition.
- ◇ Level 3 is the warmest setting and Level 1 is the least warm setting.

3.3 Heating up the product

1. Turn on the product by pressing the power button on the switch.
2. Press the on/off button to activate a different heat level.

3.4 Step-down function

The product is equipped with a step-down function:

- It switches from level 3 to level 2 after two hours.
- After another two hours, it switches from level 2 to level 1 to prevent overheating.

3.5 Charging and using the product

The product can be used both with the battery and when connected to the power outlet. Always remove the power cord from the power outlet when the product is not in use.

**NOTE**

- ◇ Actual heating time may vary depending on ambient temperature and battery condition.

4. Maintenance and cleaning

1. Turn off the product and remove the outer cover.
2. Clean the outer cover according to the washing instructions on the product washing label.
3. Clean the inner cushion with a damp sponge.
4. Do not use abrasive or harsh cleaning agents.

5. Storage

- Let the product cool down before storing it.
- Turn off the switch when not in use and store the product inside in a dry place. Store the product in a dry, clean place.
- Do not store the product in extremely high or low temperatures.
- Do not wrinkle the product by placing objects on it during storage.
- If you do not use your product for a long period of time, make sure to charge the battery at least every three months to keep the battery in good condition.

6. Disposal and recycle

Think about the environment and contribute to a cleaner living environment! You can contact your local council for all information on disposal options for discarded products.

6.1 Disposal

At the end of its life, do not dispose of the product with normal household waste, but hand it in at a government designated collection point for recycling.

Dispose of the packaging and accessories in accordance with the environmental regulations in force in your local council. Pay attention to the symbols on the different packaging materials and - if necessary - dispose of them separately. The packaging materials are marked with letters and numbers: 1-7: plastic, 20-22: paper and cardboard, 80-98 composite materials.

Separate waste into glass, paper and plastic for recycling with Tri sélectif. Reduce landfill waste and conserve natural resources. Contribute to sustainability by consciously separating.


6.2 Electric appliance

The symbol shown here with the crossed-out garbage can indicates that this product is subject to European Directive (EU) 2012/19/EU. This directive states that at the end of its service life, you must not dispose of this product with normal household waste, but hand it in at a government designated collection point for recycling.

7. Symbols

The following warnings and symbols are used in this manual, on the packaging and on the product (where applicable).

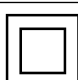Symbol	Description
	Follow the instructions in the user manual and keep it for future reference.
	A warning with this symbol indicates a potentially hazardous situation that could result in serious damage and/or fatal injury.
	A note contains additional information relevant to the use of the product.
	A warning with this symbol indicates a potentially hazardous situation that may result in minor damage and/or moderate injury.
	Wash at 30°C using the fine wash cycle.
	Do not iron.
	Do not bleach.
	Do not dry clean.
	Safety class 2 plug.

8. CE Declaration of Conformity

This product complies with all requirements and regulations according to 2014/35/EU, Directive 2014/30/EU and Directive 2009/125/EC.

The complete EU Declaration of Conformity is available at the supplier.

9. Disclaimer

Changes reserved; specifications are subject to change without stating grounds for doing so.

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©LifeGoods B.V.